


International Shipping - Live Cultures

Sample Preparation

CLIENT is responsible for sample preparation. CLIENT is responsible for culturing microbial samples and sending them to Accugenix.

Recommended shipping preparation:

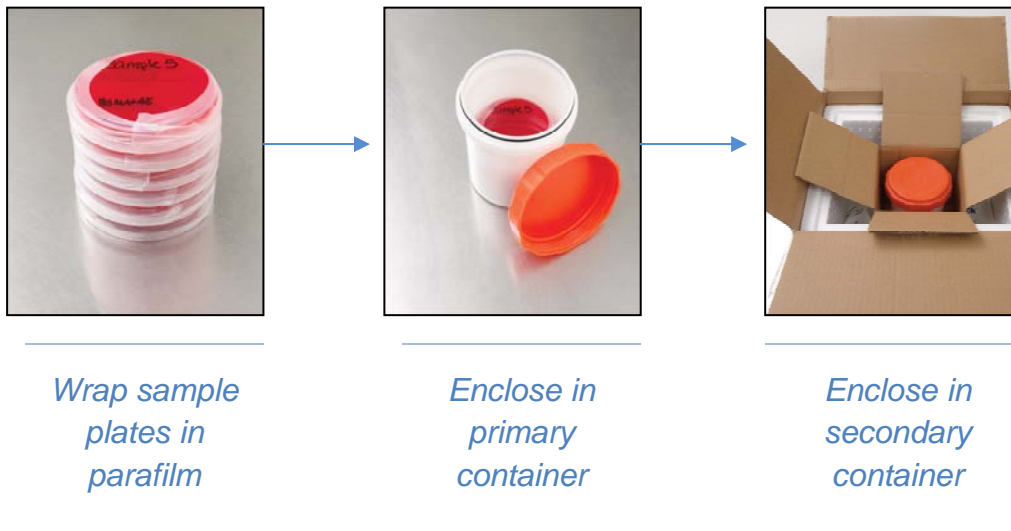

Recommendations

1. Ship at ambient temperature (packing with ice may freeze media)
2. Triple containment is recommended
 - a. Parafilm plates or slants
 - b. Protect glass from breakage by wrapping with bubble wrap
 - c. Place samples in primary container (protects secondary container if samples are damaged)
 - d. Place samples in secondary container
3. Documents for outside of package
 - a. CDC Permit
 - b. Shipper's Declaration for Dangerous Goods
 - c. Commercial Invoice or Proforma Invoice (if necessary)
 - d. International Air Waybill (FedEx, UPS, etc.)
4. Documents for inside of package
 - a. Identification Request Form for Accugenix services
 - b. Any further forms necessary for testing